




Liposuction & Fat Grafting

Infiltration Cannula

Multi-spiraled cannulas for tumescent fluid infiltration dissection prior to fat harvesting

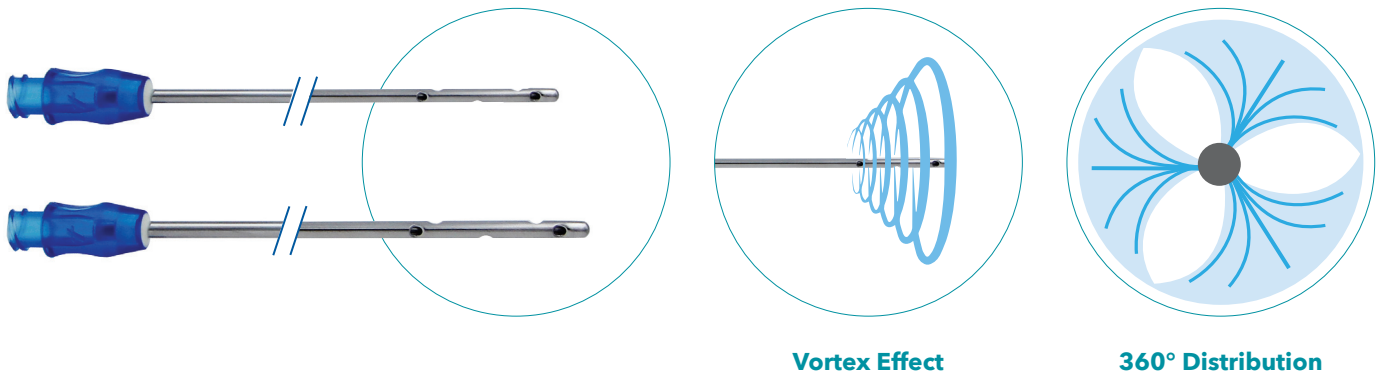









Liposuction Cannula

Body reshaping cannulas with sharp edges ports "guillotine" design to cut adipose tissue




REF	∅ mm	G	↔ mm	Holes	Box	
INFILTRATION						
	DOIN1510-2	1.5	17	100	10	10
	DOIN2015-2	2	14	150	10	10
	DOIN2525-2	2.5	13	200	10	10
	CI16150	1.65	16	150	10	10
	CI18100	1.25	18	100	10	10
	CI21100	2.1	14	100	10	10
	CI21150	2.1	14	150	10	10
	CI21250	2.5	13	250	10	10
	CI25250	2.5	13	250	10	10
	CI25300	2.5	13	300	10	10

LIPOSUCTION						
	LC31302	3	11	130	2	6
	LC32302	3	11	230	2	6
	LC42302	4	8	230	2	6
	LC42802	4	8	280	2	6
	LC31303	3	11	130	3	6
	LC32303	3	11	230	3	6
	LC42303	4	8	230	3	6
	LC42803	4	8	280	3	6
	LC52803	5	6	280	3	6
	LC62803	6	5	280	3	6
	LC32303HVS	3	11	230	3	6
	LC31304H	3	11	130	4	6
	LC32804H	4	8	280	4	6
	LC313024H	3	11	130	24	6
	LC428024H	4	8	280	24	6
	DOLF3150-2	3	11	150	3	10
	DOLF3250-2	3	11	250	3	10
	DOLF4150-2	4	8	150	3	10
	DOLF4250-2	4	8	250	3	10
	DOLF4350-2	4	8	350	3	10
	DOLF5250-2	5	6	250	3	10
	DOLF5350-2	5	6	350	3	10
	DOPL4250-2	4	8	250	7	10
	DOLP4250-2	4	8	250	8	10
	DOLP4350-2	4	8	350	8	10
	DOLF5250-2	5	6	250	8	10
	DOLP5250-2	5	6	250	8	10
	DOLP5350-2	5	6	350	8	10

Fat Harvesting Cannula

Cannulas specifically designed to perform an accurate fat SIZING during harvesting thanks to our patented diamond-shaped ports with 360° cutting edges. Macrofat Picking cannula (**REF PICK352-L**) has a luer lock hub that can also accommodate the DOTO handle to convert to liposuction via aspiration pump.

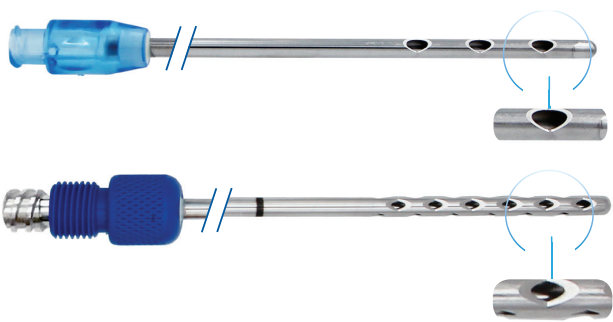

- MEGAFAT
- Small volume
- LL syringe
- Suction 360°
- MACROFAT
- Med/Large volume
- LL syringe / Pump
- Suction precision






































Metallic Luer/Lock

A range of cannula specifically designed to perform accurate fat SIZING during harvesting due to the diamond shaped cutting ports with 360° cutting edges, when selecting the appropriate negative pressure to harvest.

	REF	Ø mm	G	↔ mm	Holes	Box
HARVESTING						
	PICK3518-2	3.5	10	180	5	10
	PICK352E-2	2.0	14	150	8	10
	PICK352L-2	3.5	10	230	18	10
	PICK352A-2	4.0	8	300	8	10
	CP3015	3	11	150	2	10
	CP3023	4	8	230	2	10
	CP161053HS	1.65	16	105	3	10
	CP141253HS	2.1	14	125	3	10
	CP30153S	3	11	150	3	10
	CP141253HV	2.1	14	125	3	10
	CP31503M	3	11	150	3	10
	CP32303M	3	11	230	3	10
	CP42303M	4	8	230	3	10
	CP32303T	3	11	230	3	10
	CP42303T	4	8	230	3	10
	CP30154H	3	11	230	4	10
	CP141256HMP	2.1	14	125	6	10
	CP141506HS	2.1	14	150	6	10
	CP35176H	3	11	170	6	10
	CP141508H	3	11	150	8	10
	CP301512R	3	11	150	12	10
	CP301524H	3	11	150	24	10
	CP302324H	3	11	150	24	10
	CP122007H	2.8	12	200	7	10
	CP123007H	2.8	12	300	7	10

Fat Injection

	REF	Ø mm	G	↔ mm	Holes	Box
FAT INJECTION						
	FIL1270S-2	3.5	10	250	1	10
	FILU2515-2	2.5	12	200	1	10
	FIL2590S-2	2.7	12	200	1	10
	F9027150-2	2.7	12	150	1	10
	FIL25C5S-2	2.0	14	150	1	10
	FIL1670S-2	2.0	14	150	1	10
	F9020100-2	2.0	14	100	1	10
	F9015100-2	1.5	17	100	1	10
	F9015070-2	1.5	17	70	1	10
	F9012100-2	1.2	18	100	1	10
	F9012070-2	1.2	18	70	1	10
	CP1690V	1.65	16	90	1	10
	CP16150V	1.65	16	150	1	10
	CP1490V	2.10	14	90	1	10
	CP14150V	2.10	14	150	1	10
	MI18110KR	1.2	18	110	1	20
	MI1890KR	1.2	18	90	1	20
	MI1870KR	1.2	18	70	1	20
	MI2170KR	0.81	21	70	1	20
	MI2140KR	0.81	21	40	1	20
	MI2290KR	0.70	22	90	1	20
	MI2270KR	0.70	22	70	1	20
	MI2250KR	0.70	22	50	1	20
	MI2370KR	0.63	23	70	1	20
	MI2350KR	0.63	23	50	1	20
	MI2450KR	0.55	24	50	1	20
	MI2438KR	0.55	24	38	1	20
	MI2550KR	0.50	25	50	1	20
	MI2650KR	0.45	26	50	1	20
	MI2638KR	0.45	26	38	1	20
	MI2625KR	0.45	26	25	1	20
	MI2740KR	0.40	27	40	1	20
	POR6383	Stainless Steel Syringe Racks			18 Hole Rack	
	RAC89200	Luer Lock Transfer Adaptor			Box of 100	
	MC1711A	Luer Lock Syringe Cap			Box of 100	


st'rim®

MICROFAT grafting kit - REF PICS2130-2



SIZING



360° SUCTION

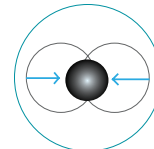

MEDIUM FLOW

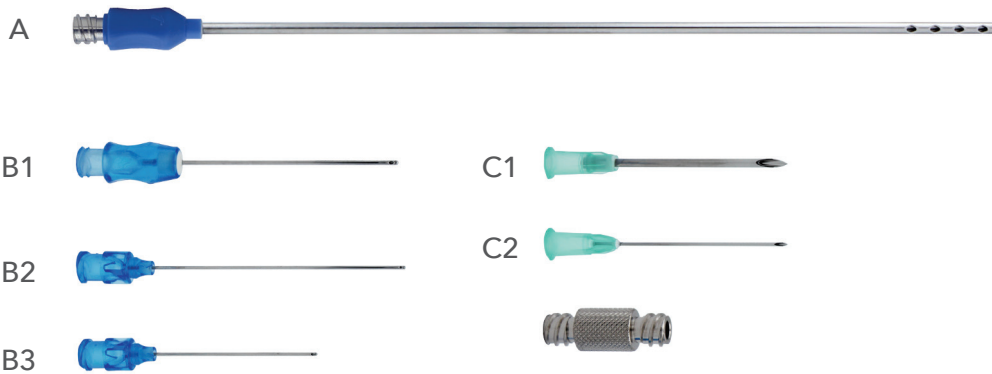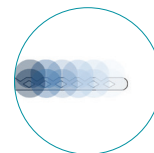

TIP	A	B1	B2	B3	C1	C2	REF	Box
∅	2.0 mm	19G	21G	21G	14G	21G	PICS2130-2	5
←mm	130	50	60	40	40	40		

	DOTO 0155-2	Detachable handle for Macrofat Picking Cannula/Macro Strim	Box of 10
--	--------------------	--	-----------

Lipocube



Nanocube "All in one" Fat Grafting - 4 ports for small volume fat grafting procedures;

Sterile, single-use device for re-sizing harvested fat particles into mini, micro and nanofat grafts for specific indications;

Supplied sterile, pack size each (1) REF L-NANO

INDICATIONS




Milifat (port 1-2) 2.4mm<

Deep fat compartments, Fat-Lifting, Bone onlay graft



Microfat (port 3) 1.2mm<

Superficial fat compartments

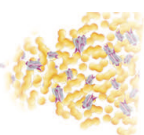

Nanofat (port 4) 600 microns <

Subdermal, Intradermal, Masotherapy, Micro-needling, topical

st'rim®

MACROFAT grafting kit - REF PICS3523-2

TIP	A	B1	B2	B3	C1	C2	REF	Box
∅	3.5 mm	12G	14G	17G	17G		PICS3523-2	5
↔mm	220	150	150	100	70			

NLF - Fat Collection System

Sterile, single use kit for fat collection during liposuction with double filter to allow swift fluid removal. There is an easy access non-return, luer lock valve for collection of high quality fat cells.

REF NLF2SYSTEM, NLF2STAND

- Completely closed and sterile circuit
- Biocompatible
- Single-use
- Compatible with the manual liposuction cannulae and power assisted devices
- Easy to use
- In compliance with the European standards - CE marking / Medical device class IIa





For post surgical support garments, visit
www.esshop.co.uk

Eurosurgical Limited
Merrow Business Centre, Guildford,
Surrey GU4 7WA, United Kingdom

+44 (0) 1483 456 007
sales@eurosurgical.co.uk
www.eurosurgical.co.uk